


# Protekt® 100

## Pressure Redistribution Foam Mattress

HCPCS  
E0184

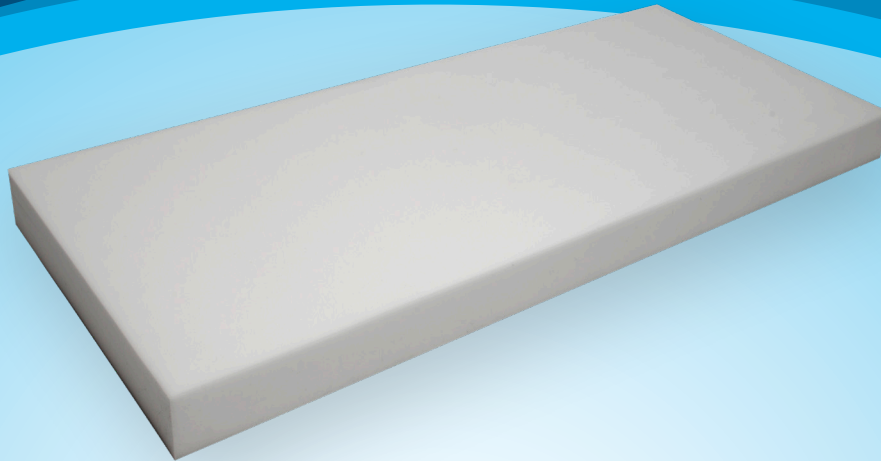

### FEATURES AND BENEFITS:

- The entire mattress core is comprised of a high resilient 1.5 density foam formulation. The ultra-soft foam construction measures at a 28 ILD (indentation load deflection) which allows for deep immersion and provides excellent pressure redistribution.
- Nylon top cover reduces perspiration and provides a low shear friction surface that is anti-microbial and fluid resistant.
- 360 degree zipper for easy removal, cleaning or replacement.
- Anti-Contamination Flap covers the zipper reducing the possibility of fluids damaging the foam core.
- Bottom Cover: Triple laminate vinyl; fluid resistant; anti-microbial/anti-bacterial; fire retardant.
- Fire sleeve meets Federal Fire Code 16 CFR part 1633.
- Weight Capacity: 350 lbs.
- Warranty: Core; 4 Year Non-Prorated. Cover; 2 Year.
- Optional raised rails available on all sizes.

81011	Protekt® 100	36"x76"x6"
81012	Protekt® 100	36"x80"x6"
81013	Protekt® 100	36"x84"x6"
81014	Protekt® 100 Plus	36"x76"x9"
81015	Protekt® 100 Plus	36"x80"x9"
81016	Protekt® 100 Plus	36"x84"x9"



270 Washington Street, Mount Vernon, NY 10553

Tel: 855-BE-PROACTIVE (855-237-7622) • Fax: 845-205-6005

www.proactivemedical.com • sales@proactivemedical.com