


Protekt[®] 300

Pressure Redistribution Foam Mattress

HCPCS
E0184

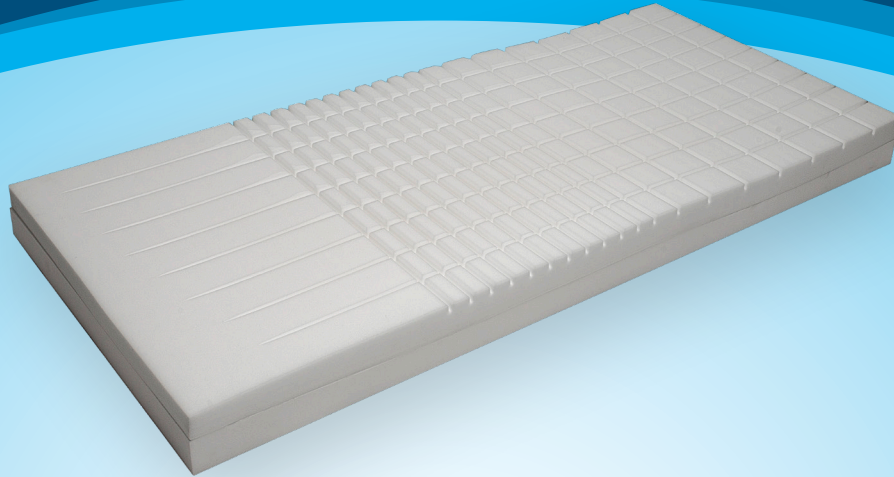

FEATURES AND BENEFITS:

- Multi-tiered pressure redistribution mattress with top layer of high density and high resilient die-cut foam in the torso and head sections offers unparalleled pressure redistribution.
- The sloped heel section combined with open cell memory foam reduces pressure and increases blood circulation by redistributing the load from heel to the less vulnerable lower leg.
- PTSS Rail System provides additional support for easy patient transfer and edge of bed sitting.
- Stretch nylon top cover reduces perspiration and provides a low shear friction surface that is anti-microbial and fluid resistant.
- 360 degree zipper for easy removal, cleaning or replacement.
- Anti-Contamination Flap covers the zipper reducing the possibility of fluid migration.
- Bottom Cover: Triple laminate vinyl; fluid resistant; anti-microbial/anti-bacterial; fire retardant.
- Fire sleeve meets Federal Fire Code 16 CFR part 1633.
- Weight Capacity: 350 lbs.
- Warranty: Core; 7 Year Non-Prorated. Cover; 2 Year.
- Optional raised rails available on all sizes.

81031	Protekt [®] 300	36"x76"x6"
81032	Protekt [®] 300	36"x80"x6"
81033	Protekt [®] 300	36"x84"x6"
81034	Protekt [®] 300 Plus	36"x76"x9"
81035	Protekt [®] 300 Plus	36"x80"x9"
81036	Protekt [®] 300 Plus	36"x84"x9"



270 Washington Street, Mount Vernon, NY 10553

Tel: 855-BE-PROACTIVE (855-237-7622) • Fax: 845-205-6005

www.proactivemedical.com • sales@proactivemedical.com